


Compression Measurement System Hohenstein
PRESSURE PROFILE
Orderer: medi Bayreuth
Item: thrombexin climax Test size: Nonstandard
Special details: size medium - thigh - adhesive welt



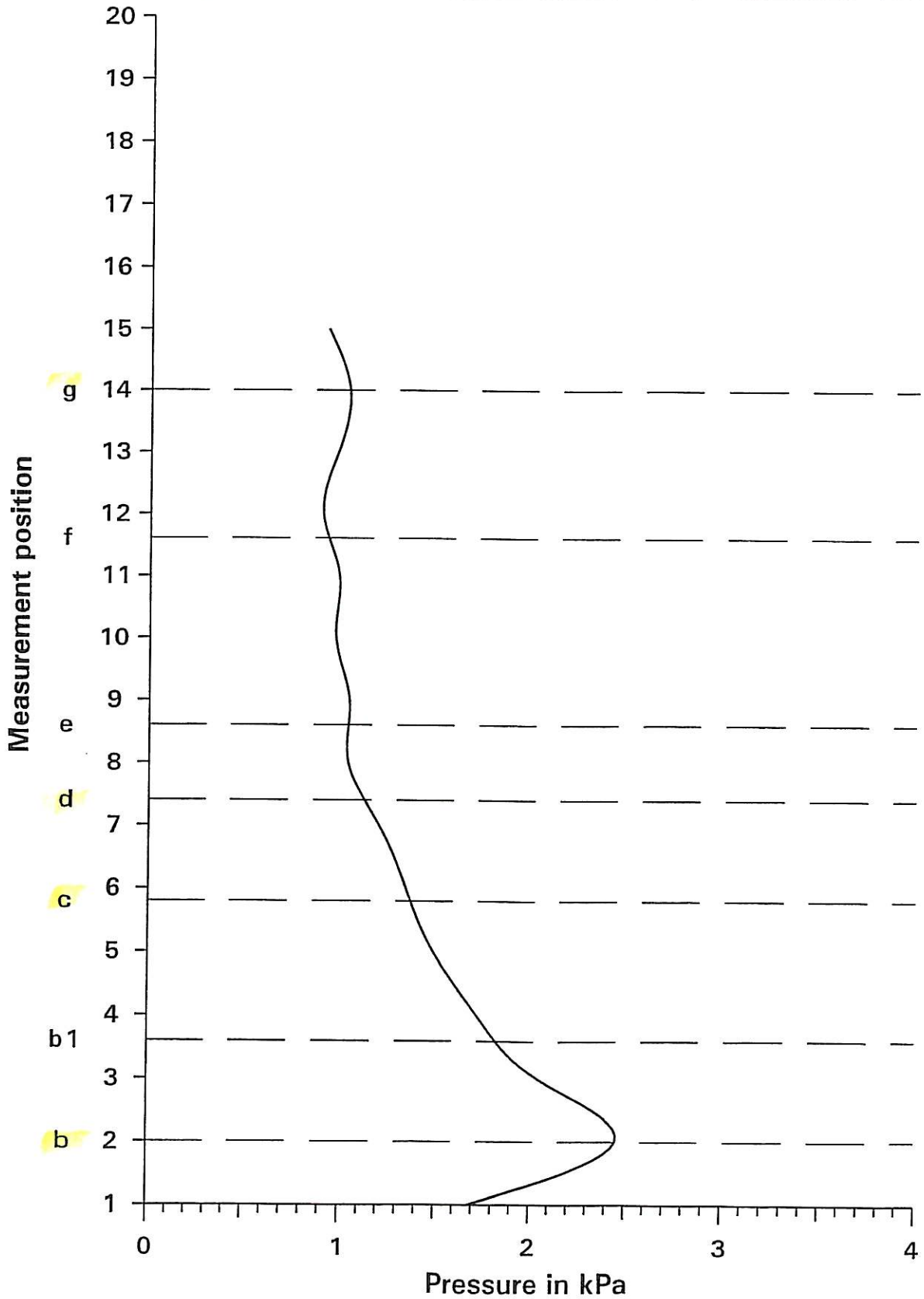

Compression Measurement System Hohenstein

TABLE OF RESULTS

Orderer: medi Bayreuth
Item: thrombexin
Manuf. size: medium
Test size: Nonstandard
Special details: thigh - with adhesive welt

Measurement Pos.	Circ. (cm)	Height (cm)	Wear stretch (%)	Tension (N/cm)	Pressure		Residual pressure (%)
					(kPa)	(mm Hg)	
b	24.0	12.0	77.6	.94	2.45	18.4	100.0
b1	30.0	20.0	85.6	.87	1.82	13.6	74.2
c	37.0	31.0	83.4	.81	1.37	10.3	55.8
d	35.0	39.0	62.6	.63	1.13	8.4	45.9
e	38.0	45.0	64.0	.63	1.04	7.8	42.4
f	55.0	60.0	97.9	.81	.93	7.0	38.1
g	63.0	72.0	66.9	1.04	1.04	7.8	42.2
1	24.6	7.0	64.0	.66	1.68	12.6	68.5
2	24.0	12.0	77.6	.94	2.45	18.4	100.0
3	27.1	17.0	86.0	.88	2.05	15.3	83.4
4	32.0	22.0	85.4	.88	1.72	12.9	70.1
5	35.9	27.0	89.2	.85	1.49	11.2	60.8
6	36.9	32.0	80.7	.79	1.34	10.1	54.8
7	35.3	37.0	65.9	.67	1.20	9.0	49.0
8	35.9	42.0	61.5	.59	1.04	7.8	42.4
9	39.9	47.0	66.6	.66	1.04	7.8	42.6
10	45.6	52.0	75.7	.71	.97	7.3	39.6
11	51.7	57.0	91.2	.81	.99	7.4	40.3
12	56.9	62.0	99.7	.82	.90	6.8	36.9
13	60.7	67.0	91.6	.95	.98	7.3	39.9
14	63.0	72.0	66.9	1.04	1.04	7.8	42.2
15	64.0	77.0	46.3	.94	.93	6.9	37.7

The results show the average of 2 measurements.